

Human IgG (Total) ELISA Kit

Cat: RK00393

This ELISA kit used for quantitation of human Total Immunoglobulin Heavy Constant Gamma (IgG (Total)) concentration in cell culture supernate, serum and plasma. For research use only, and it's highly recommended to read throughly of this manual before using the product.

Manufactured by

Global Headquarters

86 Cummings Park

Woburn, MA 01801

Tel: +8887545670

China Branch

388# Gaoxin Road (No.2)

Tel: 400-999-6126

East Lake Development Zone

E-mail: market@abclonal.com

Wuhan P. R. China

[http: www.abclonal.com.cn](http://www.abclonal.com.cn)

Contents

Introduction.....	3
Principle Of The Assay	3
Materials Provided	4
Sample Collection And Storage.....	5
Precautions For Use.....	6
Experiment Materials	7
Reagent Preparation.....	8
Wash Method.....	10
Assay Procedure	10
Assay Procedure Summary	12
Calculation Of Results.....	13
Typical Data	14
Sensitivity	15
Specificity	15
Precision	16
Recovery.....	17
Linearity Dilute	17
References.....	18

Introduction

IgG is the major immunoglobulin in blood, lymph fluid, cerebrospinal fluid, and peritoneal fluid and a key player in the humoral immune response. Serum IgG in healthy humans presents approximately 15% of total protein beside albumins, enzymes, other globulins and many more.

The Fc portion of IgG, but not F(ab')₂ or Fab fragments, can cross the placenta of a mother to enter the fetal circulation providing the fetus with postpartum protection.

IgG molecules are able to react with Fc-gamma receptors that are present on the surfaces of macrophages, neutrophils, natural killer cells, and can activate the complement system.

The binding of the Fc portion of IgG to the receptor present on a phagocyte is a critical step in the opsonizing property IgG provides to the immune response. Phagocytosis of particles coated with IgG antibodies is a vital mechanism to cope with microorganisms.

IgG is produced in a delayed response to an infection and can be retained in the body for a long time. The longevity in serum makes IgG most useful for passive immunization by transfer of this antibody. Detection of IgG usually indicates a prior infection or vaccination.

Principle Of The Assay

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for IgG (Total) has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any IgG (Total) present is bound by the immobilized antibody. Following incubation unbound samples are removed during a wash step, and then a detection antibody specific

for IgG (Total) is added to the wells and binds to the combination of capture antibody-IgG (Total) in sample. Following a wash to remove any unbound combination, and enzyme conjugate is added to the wells. Following incubation and wash steps, a substrate is added. A colored product TMB is formed in proportion to the amount of IgG (Total) present in the sample. The reaction is terminated by addition of acid and absorbance is measured. A standard curve is prepared from seven IgG (Total) standard dilutions and IgG (Total) sample concentration determined.

Materials Provided

Description	Size (192T)	Size (96T)	Size (48T)	Storage	Cat NO.
Human IgG (Total) antibody coated plate	(8×12) ×2	8×12	8×6	4°C	RM01600
Human IgG (Total) Standard lyophilized	4 vials	2 vials	1 vial	4°C	RM01597
Standard/sample Diluent (R1)	2 bottles ×20 mL	1 bottle ×20 mL	1 bottle ×20 mL	4°C	RM00023
Human IgG (Total) concentrated biotin conjugate antibody (100×)	2 vials ×120 µL	1vial ×120 µL	1 vial ×60 µL	4°C	RM01598

Biotin-Conjugate antibody Diluent (R2)	1 bottle ×32 mL	1 bottle × 16 mL	1 bottle × 16 mL	4°C	RM00024
Streptavidin-HRP concentrated (100×)	2 vials ×120 µL	1 vial ×120 µL	1 vial ×60 µL	4°C	RM01599
Streptavidin-HRP Diluent(R3)	1 bottle ×32 mL	1 bottle ×16 mL	1 bottle ×16 mL	4°C	RM00025
Wash Buffer (20×)	2 bottles ×30 mL	1 bottle × 30 mL	1 bottle ×30 mL	4°C	RM00026
Substrate Solution (Dark)	2 bottles ×12 mL	1 bottle ×12 mL	1 bottle ×6 mL	4°C	RM00027
Stop Solution	1 bottle ×24 mL	1 bottle ×12 mL	1 bottle ×12 mL	4°C	RM00028
Plate Sealers	8 strips	4 strips	2 strips		
Specification	1				

Sample Collection And Storage

1. Cell Culture Supernates:

Centrifuge 1000x g for 10 min and detect; or aliquot and store samples at -20°C to -70°C (Stored at 2-8°C if tested within 24 hours). Avoid freeze/thaw cycles.

2. Serum:

Use a serum separator tube and allow samples to clot for 30 minutes before

centrifugation for 10 minutes at 1000x g, and detect; or aliquot and store samples at -20°C to -70°C (Stored at 2-8°C if tested within 24 hours). Avoid freeze/thaw cycles.

3. Plasma

Collect plasma using EDTA or heparin as an anticoagulant. Centrifuge for 15 minutes at 1000x g within 30 minutes of collection, and detect; or aliquot and store samples at -20°C to -70°C (Stored at 2-8°C if tested within 24 hours). Avoid freeze / thaw cycles.

4. Avoid hemolytic and hyperlipidemia sample for Serum and Plasma.

5. Dilution:

Dilute samples at the appropriate multiple (recommend to do pre-test to determine the dilution factor).

Precautions For Use

- 1. Reagents may be harmful, if ingested, rinse it with an excess amount of tap water.**
- 2. Stop Solution contains strong acid. Wear eye, hand, and face protection.**
- 3. Store the kits at 2 to 8°C before use, throw away the unspent kits.**
- 4. Apart from the standard of kits, other components should not be refrigerated.**
- 5. Please perform simple centrifugation to collect the liquid before use.**
- 6. Apart from Stop Buffer and Concentrated Wash Buffer can be commonly used, the other components in the kits are specified. Do not mix or substitute reagents with those from other lots or other sources.**
- 7. Adequate mixing is very important for good result. Use a mini-vortexer at the lowest frequency.**
- 8. Mix the sample and all components in the kits adequately, and use clean**

plastainer to prepare wash buffer.

9. Both the sample and standard should be assayed in duplicate, and the sequence of the reagents should be added consistently.
10. The kit should not be used beyond the expiration date.
11. The kit should be away from light when it is stored or incubated.
12. To reduce the likelihood of blood-borne transmission of infectious agents, handle all serum, plasma and other biological fluids in accordance with NCCLS regulations.
13. To avoid cross contamination, please use disposable pipette tips.
14. Please prepare all the kit components according to the requirement. If the kits will be used several times, please seal the rest strips and preserve with desiccants. Do use up within 2 months.

Experiment Materials

1. ELIASA (measuring absorbance at 450 nm, with the correction wavelength set at 570 nm or 630 nm)
2. Pipettes and pipette tips: 0.5-10, 2-20, 20-200, 200-1000 μ L
3. Microplate washer, Squirt bottle
4. Micro-oscillator
5. Deionized or double distilled water, graduated cylinder
6. Polypropylene Test tubes for dilution

Reagent Preparation

1. Bring all reagents to room temperature before use. If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved.
2. Wash buffer: 1:20 diluted with double distilled or deionized water before use.
3. Biotin-Conjugate antibody: 1:100 diluted with the Biotin-Conjugate antibody Diluent (R2) before use, and the diluted solution should be used up within 30 min.

Dilution Method


Strip	Concentrated Biotin-Conjugate antibody (1:100)	Testing dilution buffer (R2)
2	20	1980
4	40	3960
6	60	5940
8	80	7920
10	100	9900
12	120	11880

4. Streptavidin-HRP: 1:100 diluted with the Streptavidin-HRP Diluent (R3) before use, and the diluted solution should be used up within 30 min.

Dilution Method

Strip	Concentrated Streptavidin-HRP (1:100)	Testing dilution buffer (R3)
2	20	1980
4	40	3960
6	60	5940
8	80	7920
10	100	9900
12	120	11880

5. **Standard:** Add standard/sample dilution (R1) 1mL into freeze-dried standard, sit for a minimum of 15 minutes with gentle agitation prior to making dilutions (100 ng/mL), then dilute according to the requirement (recommended concentration for standard curve: 100, 50, 25, 12.5, 6.25, 3.125, 1.5625, 0 ng/mL). Redissolved standard solution (100 ng/mL), aliquot and store at -20°C— -70°C.

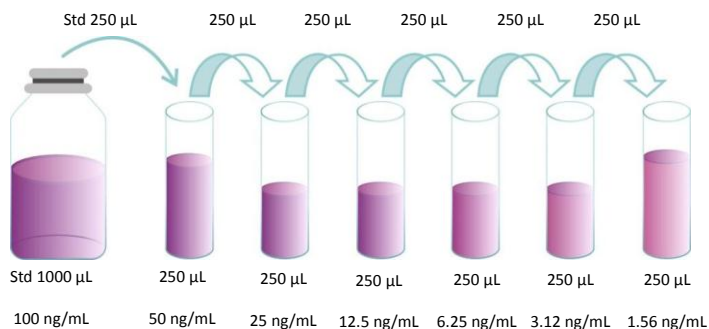

Wash Method

Automatic washer: Add wash buffer 300 μ L/well, soak for about 10-20 seconds, and wash 5 times.


Washer: Throw all the solutions in the plate well, clean with absorbent paper, and then dispense wash buffer 300 μ L/well, throw all the solutions in the plate well after holding 30 seconds, repeat 4 times.

Assay Procedure

1. Remove excess microplate strips from the plate frame, return them to the foil pouch containing the desiccant pack, and reseal.
2. Add wash buffer 300 μ L/well, aspirate each well after holding 30 seconds, repeating the process three times for a total of four washes. Then use enzyme-marked plate in a short time, do not let it dry.
3. Add 100 μ L Standard /Sample Diluent (R1) in blank well.
4. Apart from blank well, add 100 μ L different concentration of standard and sample in other wells, cover with the adhesive strip provided. Incubate for 2 hours at room temperature (20 to 25°C)
5. Wash the plate 5 times as in step 2.
6. Prepare the Biotin-Conjugate antibody Working Solution 20 minutes early.
7. Add Biotin-Conjugate antibody diluent (R2) in blank well and Biotin-Conjugate antibody Working Solution in other wells (100 μ L/well), cover with new adhesive strip provided, shake with Micro-oscillator (100 r/min). Incubate for 1 hours at room temperature (20 to 25°C)
8. Prepare the Streptavidin-HRP Working Solution 20 minutes early, place away from light at room temperature.
9. Wash the plate 5 times as in step 2.

10. Aspirate Streptavidin-HRP diluent (R3) in blank well and aspirate Streptavidin-HRP Working Solution in other wells (100 μ L/well), cover with new adhesive strip provided, shake with Micro-oscillator (100 r/min). Incubate for 30 minutes at room temperature (20 to 25°C)
11. Warm-up the ELIASA.
12. Wash the plate 5 times.
13. Aspirate substrate Solution (100 μ L/well). Incubate for 20 minutes at room temperature under dark.
14. Aspirate Stop Solution (100 μ L/well), mix, determine the optical density of each well within 5 minutes, using a microplate reader set to 450 nm. If wavelength correction is available, set to 570 nm or 630 nm. If wavelength correction is not available, subtract readings at 570 nm or 630 nm from the readings at 450nm. This subtraction will correct for optical imperfections in the plate. Readings made directly at 450 nm without correction may be higher and less accurate.

Assay Procedure Summary

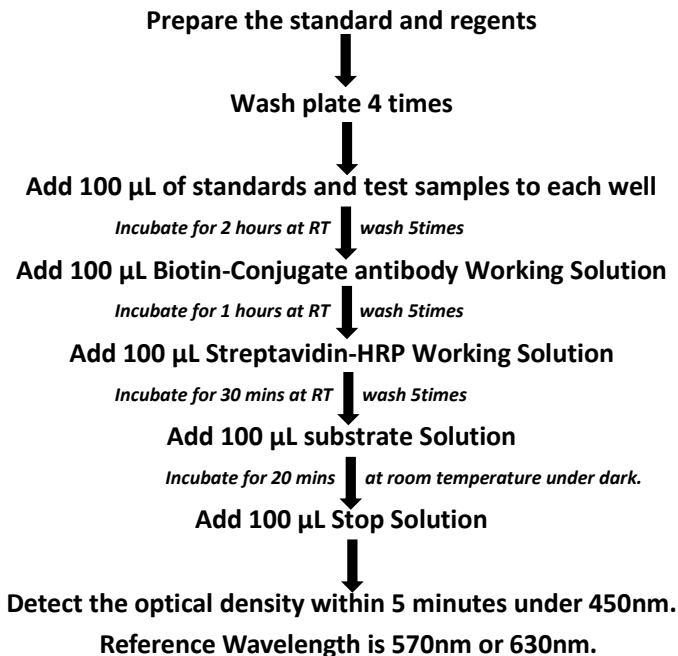


Calculation Of Results

1. Average the duplicate readings for each standard, control and sample, and subtract the average zero standard optical density (O.D.).
2. Create a standard curve by reducing the data using computer software capable of generating a log/log curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on a log/log graph. The data may be linearized by plotting the log of the IgG (Total) concentrations versus the log of the O.D. on a linear scale, and the best fit line can be determined by regression analysis.
3. If the detect result is higher than the standard curve's upper limit, then dilute samples, and the concentration read from the standard curve must be multiplied.

Typical Data

Standard (ng/mL)	OD value		Average value	Correct value
0	0.05	0.052	0.051	---
1.56	0.115	0.117	0.116	0.065
3.12	0.187	0.197	0.192	0.141
6.25	0.329	0.339	0.334	0.283
12.5	0.588	0.602	0.595	0.544
25	1.005	1.029	1.017	0.966
50	1.647	1.661	1.654	1.603
100	2.547	2.573	2.56	2.509

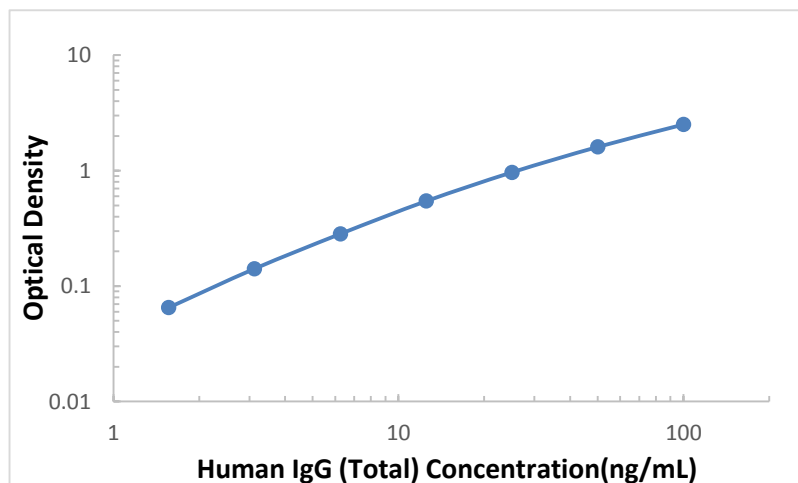

The standard curves are provided for demonstration only. A standard curve should be generated for each set of IgG (Total) assayed.

Sensitivity

The minimum detectable dose (MDD) of IgG (Total) ranged from 0.8 ng/mL. The MDD was determined by adding two standard deviations to the mean optical density value of twenty zero standard replicates and calculating the corresponding concentration.

Specificity

This assay recognizes both recombinant and natural human IgG (Total). Use 100 ng/mL to do specificity assay. No significant cross-reactivity was observed with the following:

Cross-reactive with non-human primate IgG

Precision

Intra-plate Precision

Repeat 20 times detection of 3 known concentration sample enzyme plate to evaluate the Intra-plate precision.

Sample	1	2	3
Repeat Times	20	20	20
Average Value (ng/mL)	60	310	640
Standard Deviation (SD)	3.0	13.9	26.8
Variable Coefficient CV (%)	5.0	4.5	4.2

Inter-plate Precision

Repeat 20 times detection of 3 known concentration sample enzyme plate to evaluate the Inter-plate precision.

Sample	1	2	3
Repeat Times	20	20	20
Average Value (ng/mL)	65	390	800
Standard Deviation (SD)	5.1	24.9	60
Variable Coefficient CV (%)	7.9	6.4	7.5

Recovery

Aspirate 3 different concentration of human IgG (Total) into healthy human serum and plasma, calculate the recovery.

Sample Form	Average Recover (%)	Range (%)
Serum	95	85-116
Plasma	98	83-113

Linearity Dilute

Aspirate high concentration of human IgG (Total) into 4 healthy human serum, dilute in the range of standard curve kinetics and evaluate the linearity.

Dilution	Average Value (%)	Range (%)
1:2	98	86-113
1:4	99	84-104
1:8	98	91-114
1:16	97	83-107

References

1. Fitzgerald, K., et al. Eds. 2001. The Cytokine FactsBook. Academic Press, San Diego.
2. Aggarwal, B., et al. Eds. 1992. Tumor necrosis factors: structure, function, and mechanism of action. Marcel Dekker Inc.
3. Bonavida, B., et al. Eds. 1990. Tumor necrosis factor: structure, mechanisms of action, role in disease and therapy. Karger, Basel.
4. Paul, N., et al. 1987. Annu. Rev. Immunol. 6:407.